


PAIN MANAGEMENT

PRODUCT GUIDE

neurafit™

NRFIT™ EPIDURAL RANGE

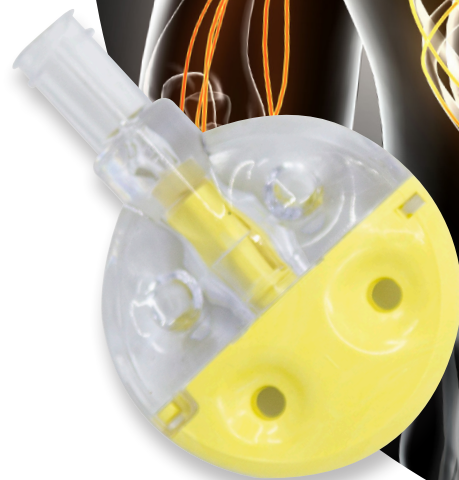

Value Life

neurafitTM

Introducing the neurafitTM NRfitTM epidural range

The new neurafitTM epidural range is NRfitTM compliant and designed with the user in mind. Easy to set up, it features the standard equipment you need to perform an epidural, with the addition of our innovative 'Twist & click connector'.

Innovative twist & click connector



Transparent section helps you see when the catheter is fully inserted.



Notches connect when closed, making it impossible to over-tighten the device and block the catheter.




Catheter is inserted fully into the device before closing to ensure a secure connection.



Audible click indicates that the device is closed and the catheter is secured.




Easy visual confirmation that the device is closed or open.

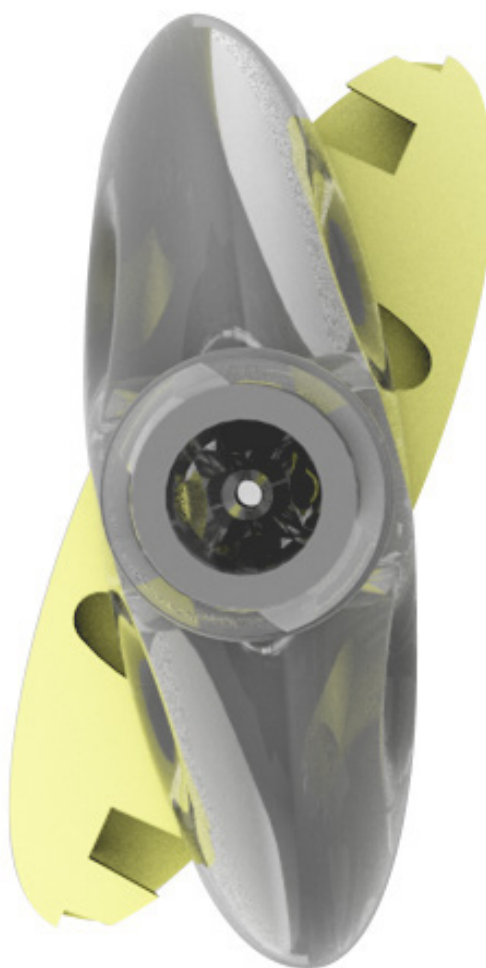

To see the twist & click connector for yourself,

scan with your smartphone

or visit vygon.co.uk/pain-airways-management/nr-fit-vygon-uk

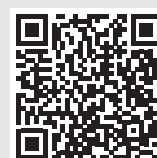

NRFit™ mini-kits

Each mini-kit includes a catheter with catheter guide, twist & click connector, Tuohy needle, 10mL LOR syringe, 0.2µm filter and yellow label.



Ordering codes

Vygon code	NHSSC Code	Description	Peripur catheter				Tuohy needle			Box/ case qty
			Size (G)	Length (mm)	Ext. ø (mm)	Int. ø (mm)	Size (G)	Length (mm)	Stylet material	
00019130A1	FTR3454	Mini-kit with peripur epidural catheter	19	90	1	0.5	17	90	Metal	20/120
00019170A1	-	Mini-kit with XRO line and peristyl catheter - a very soft transparent catheter with a metallic stylet and closed tip with 3 lateral eyes	19	90	1	0.5	17	90	Metal	20/120
0051913371A1	FTR3455	Mini-kit with XRO line and peribax catheter - a medium flex transparent catheter with no stylet	20	90	0.8	0.43	18	80	Plastic	20/120
080189322A1	FTR3453	Mini-kit with peribax catheter - a medium flex transparent catheter with a closed tip with 3 lateral eyes and no stylet	20	90	0.8	0.43	18	80	Plastic	20/120

References

- 1 Eldor, J., To NRFit or not to NRFit?, J Health Sci Dev Vol: 2, Issu: 1 (09-13), 2019
- 2 Yentis, S. M., Randall K., Drug errors in obstetric anaesthesia: a national survey, International Journal of Obstetric Anaesthesia, Volume 12, Issue 4, 246 - 249
- 3 ISO 80369-1:2018 (en) Small-bore connectors for liquids and gases in healthcare applications — Part 1: General, requirements <https://www.iso.org/obp/ui/#iso:std:iso:80369:-1:ed-2:v1:en>, visited on 16/12/2019
- 4 ESRA Congress 2019 Questions form.

09/2023 - DXJB0035008 V00657 v6

FOR FURTHER INFORMATION, PLEASE CONTACT: info@vygon.co.uk

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature. This brochure has been printed on responsibly sourced and sustainable material. To help us reduce our carbon footprint all of our literature is available electronically either from your Product Specialist, on our website or by emailing info@vygon.co.uk

VYGON (UK) LTD, THE PIERRE SIMONET BUILDING, V PARK,
GATEWAY NORTH, LATHAM ROAD, SWINDON, WILTSHIRE SN25 4DL
RECEPTION: +44 (0)1793 748800 WWW.VYGON.CO.UK

    [@vygonuk](https://www.instagram.com/vygonuk)


Value Life